




→ PIEZOSURGERY® MEDICAL
SURGICAL INSERTS

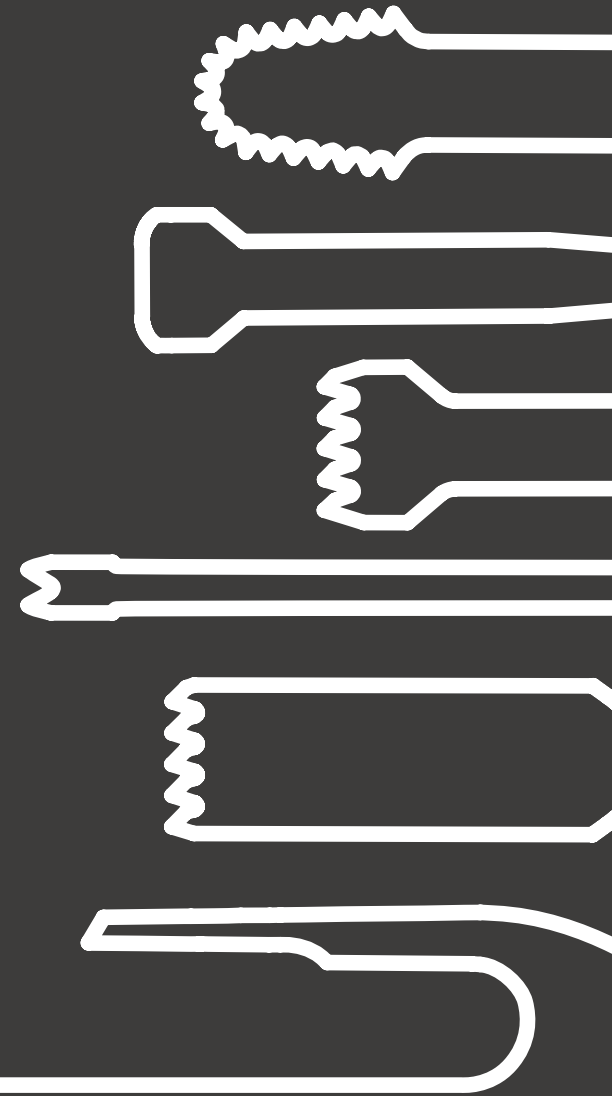

SURGICAL INSERTS

→ EXCELLENCE IN QUALITY



Intraoperative advantages and postoperative benefits are dependant on the quality of the osteotomies performed. Passion, dedication, and expertise are at the forefront of our insert production process.

We use only medical grade stainless steel in the production of Mectron inserts. Every single ultrasonic insert has to pass 12 working steps before it is ready to bear our name.



watch the inserts
development process!

→ HOW ARE THE INSERTS DEVELOPED AND PRODUCED?

Our inserts are made in Italy, in our
production site based in the Liguria Region.

SURGICAL INSERTS

→ EXCELLENCE IN VARIETY

Osteotomy, Osteoplasty, Drilling, Finishing – PIEZOSURGERY® medical inserts cover a vast variety of surgical needs. Whatever your choice is, there is one thing they all have in common: they offer the best performance you will find on the market.

Based on the specific surgical application, each insert is designed to work with a specific handpiece: standard or plus, that differ in terms of power and ergonomics.

Surgical inserts are SINGLE USE. This way each insert is always sharp and effective, ensuring each osteotomy is precise and accurate.



SURGICAL INSERTS

→ FOR STANDARD HANDPIECE

The standard handpiece has been designed to be used in osteotomies where precision and cutting control make the difference.

The “gray” handpiece features the widest range of inserts for all the four different bone cutting applications.



→ OSTEOTOMY

Surgical inserts of different shapes and dimensions, short and long, curved and angled, designed to perform osteotomies with the utmost safety even in difficult to reach surgical sites.

- Saw thickness from 0.35 to 0.6 mm
- Osteotomy depth up 20 mm



→ OSTEOPLASTY

Surgical inserts curved and angled, with sharp edges, for bone modeling and bone chip harvesting.



→ DRILLING

Surgical inserts to drill holes with very tight tolerance, minimizing the risk of bone necrosis.

- Head diameters from 0.8 to 1.8 mm



→ FINISHING

Surgical inserts of different shapes and dimensions, curved and angled, with heads of different shapes and with different diamond coatings, to finish the osteotomies in very delicate anatomies.



→ OSTEOTOMY SURGICAL INSERTS



→ MT1-10

Osteotomy micro saw
Operative length: 10.7 mm
Saw width: 3.8 mm
Saw thickness: 0.55 mm

→ MT1S-10

Osteotomy micro saw
Operative length: 10.7 mm
Saw width: 3 mm
Saw thickness: 0.35 mm

→ MT1-20

Osteotomy micro saw
Operative length: 20 mm
Saw width: 3.8 mm
Saw thickness: 0.6 mm

→ MT9-13

Osteotomy micro saw
Length insert: 42.5 mm
Operative length: 12 mm
Saw width: 3.3 - 4.4 mm
Saw thickness: 0.35 mm

→ UNIVR

**Round shape
osteotomy micro saw**
Length insert: 41.9 mm
Operative length: 5 mm
Saw width: 4.5 mm
Saw thickness: 0.5 mm

→ MT6S-10

Osteotomy micro saw
Operative length: 10 mm
Saw width: 3.6 mm
Saw thickness: 0.35 mm

→ OSTEOTOMY SURGICAL INSERTS



→ MT2L-4

Left angled micro saw
Operative length: 4 mm
Saw width: 3.8 mm
Saw thickness: 0.6 mm

→ MT2R-4

Right angled micro saw
Operative length: 4 mm
Saw width: 3.8 mm
Saw thickness: 0.6 mm

→ MT1204

Angled micro saw
Operative length: 3.5 mm
Saw width: 4.2 mm
Saw thickness: 0.55 mm

→ MT11L-10

Left angled micro saw
Operative length: 10 mm
Saw width: 4 mm
Saw thickness: 0.55 mm

→ MT11R-10

Right angled micro saw
Operative length: 10 mm
Saw width: 4 mm
Saw thickness: 0.55 mm

→ MT7-3

Osteotomy micro saw
Length insert: 44.6 mm
Operative length: 3 mm
Saw width: 3.6 mm
Saw thickness: 0.8 mm

→ OSTEOPLASTY SURGICAL INSERTS



→ MP1

Osteoplasty trapezoidal scalpel

Length insert: 27.6 mm
Operative length: 2.4 mm
Saw width: 2.8>3.7 mm
Saw thickness: 0.8 mm



→ MP2

Osteoplasty circular scalpel

Length insert: 27.3 mm
Operative length: 4 mm
Saw width: 3.9 mm
Saw thickness: 0.8 mm



→ MP3-a30

Osteoplasty circular scalpel

Length insert: 35.6 mm
Operative length: 3 mm
Saw width: 2.4 mm
Saw thickness: 0.8 mm
Operative angle: 29°



→ MP0745

Sharped semi sphere

Operative length: 1.1 mm
Saw width: Ø 4.5 mm

→ DRILLING SURGICAL INSERTS



→ MD2-08

Micro-perforation
sharp cylinder insert
Operative length: 12 mm
Ø 0.8 mm



→ MD2-10

Micro-perforation
sharp cylinder insert
Operative length: 12 mm
Ø 1.0 mm



→ MD3-14

Micro-perforation
sharp cylinder insert
Operative length: 12 mm
Ø 1.4 mm



→ MD3-18

Micro-perforation
sharp cylinder insert
Operative length: 12 mm
Ø 1.8 mm

→ FINISHING SURGICAL INSERTS



→ MF1

Diamond coated flap scalpel
Operative length: 3 mm
Operative width: 3.7 mm
Operative thickness: 0.85 mm
Diamond grit size (µm): D91

→ MF2

Diamond cylinder insert
Operative length: 5.3 mm
Operative width: Ø 2.4 mm
Diamond grit size (µm): D151

→ MF3

Diamond ball insert
Operative length: 1.7 mm
Operative width: Ø 1.7 mm
Diamond grit size (µm): D91

→ MF4

Diamond conic insert
Operative length: 10 mm
Operative width: Ø 0.6>1.5 mm
Diamond grit size (µm): D30

→ MF5

Diamond truncated cone insert
Operative length: 13 mm
Operative width: 2.1>2.7 mm
Operative thickness: 1.1>1.9 mm
Diamond grit size (µm): D126

→ MF6

Blunt cone compressor
Operative length: 0.5 mm
Operative width: 4.2 mm

SURGICAL INSERTS

→ FOR PLUS HANDPIECE


The plus handpiece has been designed to be used in osteotomies of highly mineralized bone, where more power and precision are required.

The “black” handpiece features different type of inserts for all the three different bone cutting applications. Long inserts with an extension and a cover to preserve delicate structures and soft tissues are available.

→ OSTEOTOMY

Surgical inserts of different shapes and dimensions, short and long, curved and angled, designed to perform osteotomies with the utmost safety even in difficult to reach surgical sites.

- Saw thickness from 0.35 to 0.6 mm
- Osteotomy depth up to 20 mm
- Shank length up to 10 cm

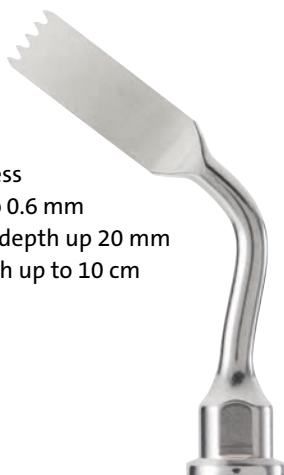

→ OSTEOPLASTY

Surgical inserts short and long, curved and angled, with sharp edges, for bone modeling and bone chip harvesting.

- Shank length up to 10 cm



→ FINISHING

Surgical inserts with diamond coating, to finish the osteotomies in very delicate anatomies.



→ OSTEOTOMY SURGICAL INSERTS



→ MT4-10 +

Osteotomy micro saw
Operative length: 12 mm
Saw width: 4.1 mm
Saw thickness: 0.55 mm

→ MT4-20 +

Osteotomy micro saw
Operative length: 19.3 mm
Saw width: 4.2 mm
Saw thickness: 0.55 mm

→ MT10-20 +

Osteotomy micro saw
Operative length: 19.8 mm
Saw width: 4.2 mm
Saw thickness: 0.55 mm

→ MT1412 +

Osteotomy micro saw
Operative length: 12 mm
Saw width: 3.3>4.4 mm
Saw thickness: 0.35 mm

→ MT1520 +

Osteotomy micro saw
Operative length: 20 mm
Operative width: 2.95>4.8 mm
Operative thickness: 0.45 mm

→ OSTEOTOMY SURGICAL INSERTS



→ MT8-20 L

Long osteotomy micro saw
 Operative length: 19.1 mm
 Saw width: 2.95>4.26 mm
 Saw thickness: 0.45 mm
 Total length insert + extension:
 117 mm



→ MT5-10 L

Long osteotomy micro saw
 Operative length: 13 mm
 Saw width: 3 mm
 Saw thickness: 0.55 mm
 Total length insert + extension:
 116.3 mm



→ OSTEOPLASTY SURGICAL INSERTS



→ MP4 +

Osteoplasty circular scalpel
Operative length: 4 mm
Saw width: 3.9 mm
Saw thickness: 0.8 mm



→ MP5 L

Long osteoplasty scalpel
Operative length: 6.8 mm
Saw width: 3.9 mm
Saw thickness: 0.8 mm
Total length insert + extension: 112.8 mm



→ MP6 L

Long and straight osteoplasty scalpel
Operative length: 6.4 mm
Saw width: 3 mm
Saw thickness: 0.9 mm
Total length insert + extension: 114.8 mm



→ FINISHING SURGICAL INSERTS



→ MF0702 +

Diamond coated ball insert
Operative length: 10.2 mm
Ø 1.7 mm
Diamond grit size (µm): D90



→ MF0804 +

Diamond coated ball insert
Operative length: 3.6 mm
Ø 3.6 mm
Diamond grit size (µm): D150



→ MF1004 L

Diamond coated ball insert
Operative length: 43.4 mm
Ø 3.5 mm
Total length insert + extension:
120 mm



PIEZOSURGERY® IN PRACTICE

mectron
medical technology



YOUTUBE



FACEBOOK






LINKEDIN



INSTAGRAM

Would you like to learn more about clinical protocols and the experiences of surgeons who collaborate with us? Discover the world of PIEZOSURGERY® and join the community on our social channels and website. You will have at your disposal:

- CLINICAL ANIMATIONS
- SURGICAL PROTOCOLS
- WEBINARS
- USERS TESTIMONIAL

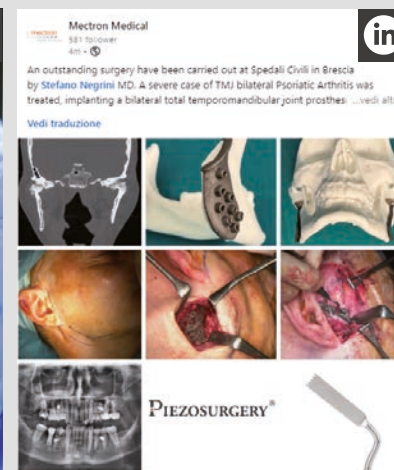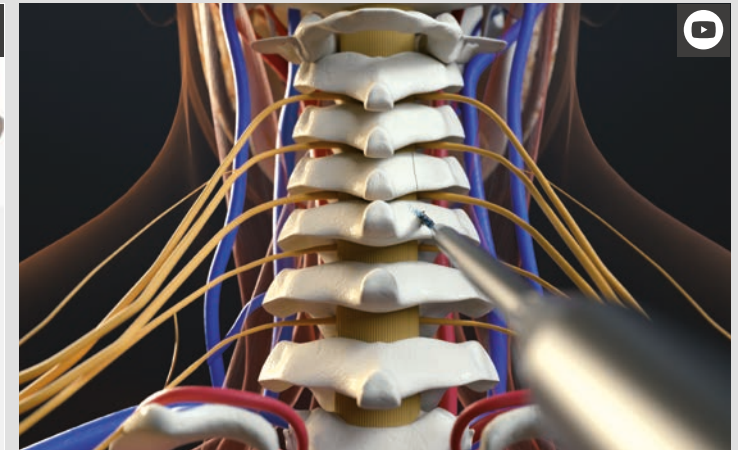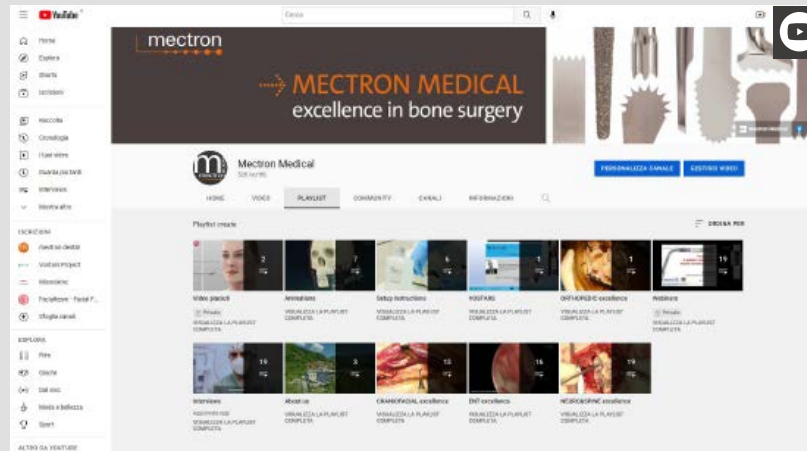

PIEZOSURGERY® DEVICES

→ EXCELLENCE IN BONE SURGERY

PIEZOSURGERY® *plus* is the complete device: it is engineered to support an expansive range of surgical applications.

→ COMPLETE

- Power joins precision in two handpieces
- High efficiency
- Cutting-edge technology



→ SURGICAL APPLICATIONS

ORAL/MAXILLOFACIAL SURGERY
OTOLARYNGOLOGY
PLASTIC/RECONSTRUCTIVE
HAND SURGERY
FOOT SURGERY
NEUROSURGERY
SPINE SURGERY
ORTHOPEDICS
THORACIC SURGERY

→ PIEZOSURGERY® *plus*

05170003

→ PIEZOSURGERY® plus

→ PIEZOSURGERY® flex



PIEZOSURGERY® *flex* is the essential device: with its perfect control, precision and ease of use it is your perfect start into a new era of bone surgery.

→ ESSENTIAL

- Control joins precision in one handpiece
- Compact and portable
- Easy to use



→ PIEZOSURGERY® *flex*

05170002



PRODUCTS

ACCESSORIES

Handpiece for plus channel (WAVE)	03120249
Torque wrench for plus channel	02900116
Torque wrench for extension of long inserts	02900115
Handpiece for standard channel	03120127
Torque wrench for standard channel	02900080
Irrigation kit single use (box of 10 units)	03230008
Trolley-case	04440018
Cart	03540009

CONTAINERS FOR STERILIZATION OF HANDPIECE AND TORQUE WRENCH

Container for sterilization completed with cover with valve	02900173
Tray for sterilization completed with cover	02900172
Thermodisinfection adaptor for handpiece	04610008
Filter for thermodisinfection adaptor	04590006

SPARE PARTS

Power-supply cable	00050020
Footswitch for PS <i>plus</i>	04620004
Peristaltic pump	03210006
Drip stands for saline bag	01380002
Metal cone for standard handpiece	00710034
Metal cone for plus handpiece	00710090
Protection for standard handpiece's connector	03150086
Protection for plus handpiece's connector	00690022

→ SURGICAL INSERTS

→ STANDARD HANDPIECE 03120127

→ PLUS HANDPIECE 03120249

OSTEOTOMY

MT1-10	03600001
MT1S-10	03600007
MT1-20	03600002
MT1204	03600019
MT2R-4	03600003
MT2L-4	03600004
UNIVR	03600008
MT6S-10	03600011
MT7-3	03600012
MT9-13	03600016

MT4-10 +	03600010
MT5-10 L	03600009
MT8-20 L	03600013
MT4-20 +	03600014
MT10-20 +	03600015
MT1412 +	03600020
MT1520 +	03600021

OSTEOPLASTY

MP1	03610001
MP2	03610002
MP3-a30	03610003
MP0745	03610010

MP4 +	03610007
MP5 L	03610008
MP6 L	03610009

DRILLING

MD2-08	03620010
MD2-10	03620004
MD3-14	03620006
MD3-18	03620008

FINISHING

MF1	03630001
MF2	03630002
MF3	03630003
MF4	03630004
MF5	03630005
MF6	03630006

MF0702 +	03630007
MF0804 +	03630008
MF1004 L	03630009



mectron s.p.a.

via Loreto 15/A, 16042 Carasco (Ge), Italy

tel +39 0185 35361, fax +39 0185 351374

→ www.mectron.com – piezosurgery@mectron.com



© Copyright mectron S.p.A., Carasco, Italy

All rights reserved. Texts, pictures and graphics of mectron brochures are protected by copyright and other protection laws.

Without written approval of mectron S.p.A. the contents may not be copied, distributed, changed or made available to third parties for commercial purposes.

The images shown are for illustration purposes only and may not be an exact representation of the product. Please consult product labels and manual for indications, contraindications, hazards, warnings, cautions.

Please contact your local sales representative for product availability.

All products illustrated in this brochure are CE-marked and registered in accordance with the EU Medical Device Regulation (MDR 2017/745).

These medical devices are intended exclusively for use by qualified healthcare professionals. The content of this brochure is intended for healthcare professionals for informational and educational purposes only and does not constitute medical advice. Please refer to the Instructions for Use (IFU) and local regulations prior to use.

Get in touch with our team!

